

Your practice is *our* inspiration.™

**Maxcem Elite** | Simplicity refined



# Maxcem Elite

Taking simplicity to the next level.



Maxcem Elite self-etch, self-adhesive resin cement for indirect restorations is simplicity refined—with higher bond strengths and superior ease of use—offering proven technology without sacrificing performance. Easy cleanup, no hand mixing, and no refrigeration required. For more convenience, intraoral and root-canal tips allow for dispensing of material into hard-to-reach places.

*Inspired by:  
Bonding agents invite  
technique sensitivity, a  
daily challenge for dentists.  
Maxcem Elite eliminates this  
step for great performance,  
convenience—and peace  
of mind.*

## Research

Maxcem Elite™ employs the same unique adhesive used in Kerr's OptiBond family of adhesives. The optimized resin matrix and filler systems improve wettability, resulting in greater adhesion and allowing for immediate light curing and better shelf-life stability. Maxcem Elite's thixotropic property facilitates the placement of material and easy cleanup. A proprietary redox initiator system eliminates the inherent discoloration of BPO/tertiary amine initiator systems for a more esthetic restoration.


## Features

**Bond strengths between 22-36 MPa.** Depending on substrates. High bond strengths without an adhesive.


**Ease of use.** No need to refrigerate, easy cleanup, and no hand mixing required.

**Compatibility with common substrates.** Anterior or posterior ceramics, PFMs and metal restorations, and CAD/CAM materials.

**Efficient self-cure mechanism.** Sets up completely in the absence of light.

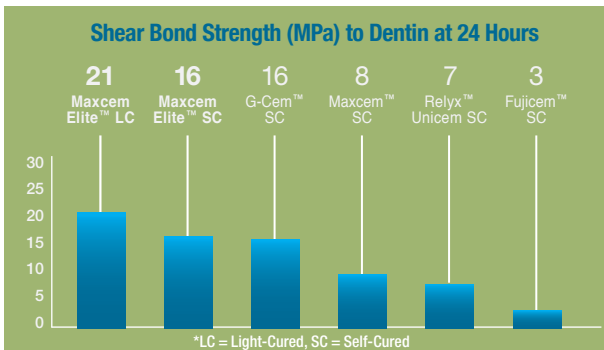

24-hour testing. Data from Dr. Carlos Muñoz, SUNY(Buffalo). Available upon request.

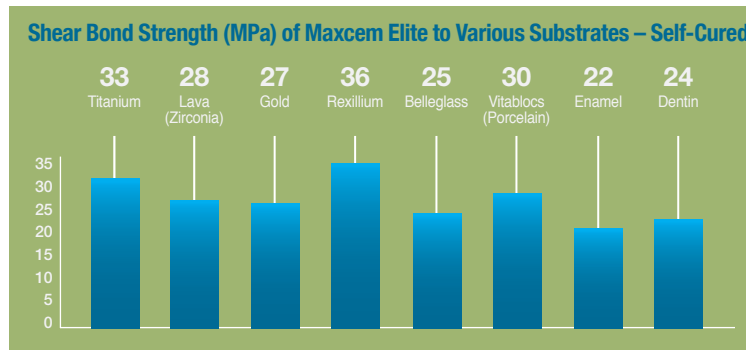

24-hour testing. Internal data. Available upon request.

## Ordering

### 34417 Maxcem Elite Intro Kit

Contains: 3 dual syringes (5 g): 2 clear, 1 white; 20 4:1 automix tips (regular); 10 4:1 automix tips (wide); 10 intraoral tips; 10 root-canal tips; Directions for use; Technique guide

### Maxcem Elite Individual Shade Refill Packs 5 Shades

Contains: 2 dual syringes (5 g) choice of one shade; 10 4:1 automix tips (regular); 10 4:1 automix tips (wide); 10 intraoral tips; Directions for use; Technique guide

- 33872 Clear
- 33873 White
- 33874 White Opaque
- 33875 Yellow
- 33876 Brown

### 34418 Maxcem Elite Bulk Pack

Contains: 4 dual syringes Maxcem Elite (5 g): clear; 20 4:1 automix tips (regular); 10 4:1 automix tips (wide); 10 intraoral tips; 10 root-canal tips; Directions for use

### Accessories

- 32862 50 4:1 automix tips (regular)
- 33879 10 4:1 automix tips (wide)
- 33880 10 intraoral tips
- 33881 10 root-canal tips



©2009 Kerr Corporation Lit. No. 34081 Rev. 2  
All trademarks are property of their respective owners.